

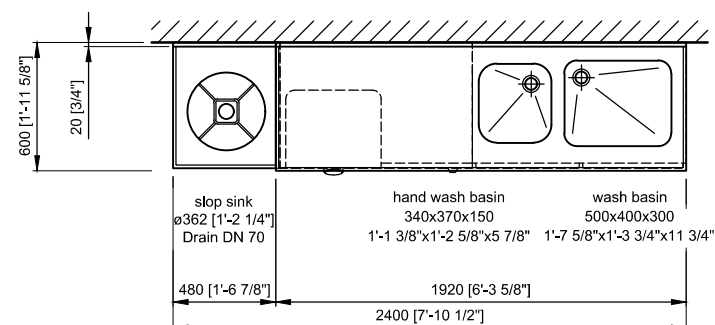




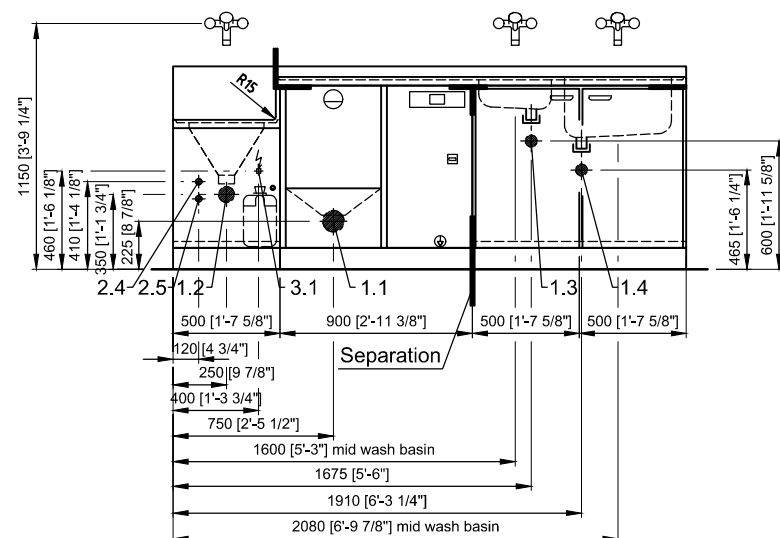


Wall connection detail X

Table top with broad upstand, straight with hygiene radius
S 1:2

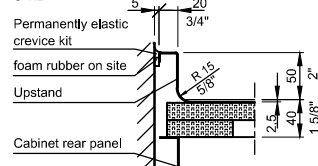

Slop sink cover

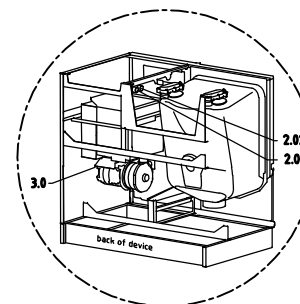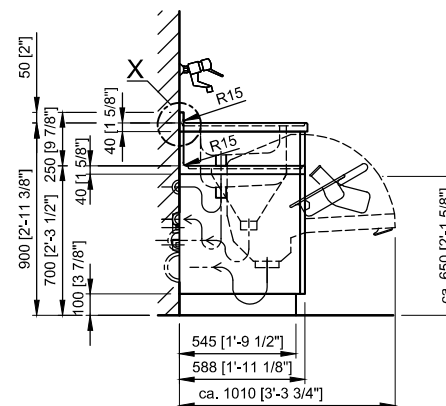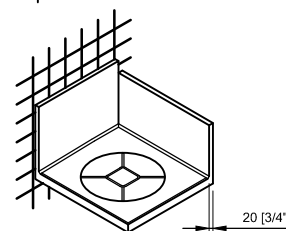

Appliance connections

MEIKO TopLine

1.0	Drain on the machine: PE siphon, $\varnothing$ 110 mm
1.1	On-site drain: Wall sleeve, DN 100
1.2	Drain on slop sink unit: PE odour trap, $\varnothing$ 75 mm On-site drain: Wall sleeve, DN 70
1.3	Drain on hand wash basin G 1 1/2" Drain on site: odour trap $\varnothing$ 40/ $\varnothing$ 50 mm, wall sleeve DN50 On site: drill hole in cabinet rear wall
1.4	Drain on wash basin G 1 1/2" Drain on site: odour trap $\varnothing$ 40/ $\varnothing$ 50 mm, wall sleeve DN50 On site: drill hole in cabinet rear wall
2.01	Hot water connection on the machine: G 3/4"a, dirt trap
2.4	On-site water connection: hot water, 45-65°C with stop valve Flow: 18 l/min, flow pressure: 100-500 kPa
2.02	Cold water connection on the machine: G 3/4"a, dirt trap
2.5	On-site water connection: cold water, 5-25°C with stop valve Flow: 18 l/min, flow pressure: 100-500 kPa
3.0	Electrical connection on unit: 3 Gnd ~ 208/230V, 60Hz, 3,9/4,6kW ⊕ Potential equalisation (M6 bolt) Protection class: Equipment Class 1, Pollution Degree 2, Installation Category II. Main switch. Connect cable approx. 2,5m
3.1	On-site electrical connection: Local fuse protection: 3 x 15A Protective measures and connection of the equipotential bonding must be carried out according to IEC 364-5-54 and conform to the local and national regulations. If the appliance is hard wired, a multiple-pole main switch must be installed. The main switch must be compliant with IEC 60947-1 and IEC 60947-3. Install the main switch in the vicinity of the appliance and in an area that is easily accessible. It must be clear that the main switch is linked to the appliance.

attachment "important information" to be considered!

Machine equipment

DT-ND/V Decentralised thermal disinfection

Slop sink flushing via MEIKO TopLine

-	Original	31.07.2025	rgd
	MEIKO USA, INC. 1349 HEIL QUAKER BLVD, LA VERGNE, TN 37086 UNITED STATES OF AMERICA TEL: (615) 399-6600 FAX: (615) 399-6620 WEB: www.meiko.us © 2017 MEIKO USA, Inc. All rights reserved.	It is neither allowed to give this drawing without our permission to a third party for information or reproduction, nor may it be used for competition purposes. We reserve all rights !	
Reference: - -		Type: SAN24BW/MEIKO TOPLINE40 M2 DT-ND/V	
Drawing-No.: P00052298		Project-No.: - KD M2 USA	Order-No. KG: Order-No.:
Scale: 1:25	drawn: 31.07.25 rgd	checked: oq-mt	Customer checked: